

3M™ Surgical Tapes & Dressings

3M™ Surgical Tapes & Dressings

Cover, care and protect

Ordering Information

Product	Ref no.	Dimensions	No. per box	
3M Cloth Adhesive Tape	2950-0	1,25 cm x 9,14 m	24 rolls	
	2950-1	2,5 cm x 9,14 m	12 rolls	
	2950-2	5,0 cm x 9,14 m	6 rolls	
	2950-3	7,5 cm x 9,14 m	4 rolls	
	Also available in 5 metre length			
	2952M	3,8 cm x 10 m	12 rolls	
With Spool	2951S-0	1,25 cm x 5 m	32 rolls	
	2951S-1	2,0 cm x 5 m	28 rolls	
	2951S-2	5,0 cm x 5 m	12 rolls	
3M Durapore™ Tape	1538-0	1,25 cm x 9,14 m	24 rolls	
	1538-1	2,5 cm x 9,14 m	12 rolls	
	1538-2	5,0 cm x 9,14 m	6 rolls	
	1538-3	7,5 cm x 9,14 m	4 rolls	
3M Transpore™ White Tape	1534-0	1,25 cm x 9,14 m	24 rolls	
	1534-1	2,5 cm x 9,14 m	12 rolls	
	1534-2	5,0 cm x 9,14 m	6 rolls	
	1534-3	7,5 cm x 9,14 m	4 rolls	
	Also available in 5 metre length			
3M Transpore™ Tape	1527-0	1,25 cm x 9,14 m	24 rolls	
	1527-1	2,5 cm x 9,14 m	12 rolls	
	1527-2	5,0 cm x 9,14 m	6 rolls	
	1527-3	7,5 cm x 9,14 m	4 rolls	
	Also available in 5 metre length			
3M Micropore™ Tape	1530-0	1,25 cm x 9,14 m	24 rolls	
	1530-1	2,5 cm x 9,14 m	12 rolls	
	1530-2	5,0 cm x 9,14 m	6 rolls	
	1530-3	7,5 cm x 9,14 m	4 rolls	
	Also available in 5 metre length			
	With dispenser	1535-0	1,25 cm x 9,14 m	24 rolls
		1535-1	2,5 cm x 9,14 m	12 rolls
		1535-2	5,0 cm x 9,14 m	6 rolls
		1535-3	7,5 cm x 9,14 m	4 rolls
	Tan	1533-0	1,25 cm x 9,14 m	24 rolls
1533-1		2,5 cm x 9,14 m	12 rolls	
1533-2		5,0 cm x 9,14 m	6 rolls	
3M Blenderm™ Tape	1525-0	1,25 cm x 4,6 m	24 rolls	
	1525-1	2,5 cm x 4,6 m	12 rolls	
	1525-2	5,0 cm x 4,6 m	6 rolls	
3M Microfoam™ Tape	1528-1	2,5 cm x 5,0 m	12 rolls	
	1528-2	5,0 cm x 5,0 m	6 rolls	
	1528-3	7,5 cm x 5,0 m	4 rolls	
	1528-4	10,0 cm x 5,0 m	3 rolls	

Product	Ref no.	Dimensions	No. per box	
3M Medipore™ H Perforated Tape	2861	2,5 cm x 9,14 m	24 rolls	
	2862	5,8 cm x 9,14 m	12 rolls	
	2863	7,5 cm x 9,14 m	12 rolls	
	2864	10,0 cm x 9,14 m	12 rolls	
	2866	15,0 cm x 9,14 m	12 rolls	
	2868	20,0 cm x 9,14 m	6 rolls	
3M Medipore™ Tape with liner	2991-1	5 cm x 10 m	1 roll	
	2991-2	10 cm x 10 m	1 roll	
	2991-3	15 cm x 10 m	1 roll	
	2991-4	20 cm x 20 m	1 roll	
	2991-5	30 cm x 10 m	1 roll	
	2991-6	2,5 cm x 10 m	1 roll	
3M Medipore™ + Pad Dressing	3562	5 cm x 7,2 cm	50 dressings	
	3564	6 cm x 10 cm	50 dressings	
	3566	10 cm x 10 cm	25 dressings	
	3569	10 cm x 15 cm	25 dressings	
	3570	10 cm x 20 cm	25 dressings	
	3571	10 cm x 25 cm	25 dressings	
	3573	10 cm x 35 cm	25 dressings	
3M Tegaderm™ + Pad Dressing with absorbent pad	3582	5 cm x 7,2 cm	50 dressings	
	3584	6 cm x 10 cm	50 dressings	
	3586	9 cm x 10 cm	25 dressings	
	3587	9 cm x 10,5 cm	25 dressings	
	3589	9 cm x 15 cm	25 dressings	
	3590	9 cm x 20 cm	25 dressings	
	3591	9 cm x 25 cm	25 dressings	
	3593	9 cm x 35 cm	25 dressings	
	Pharmacy packs	3582P	5 cm x 7,2 cm	5 dressings
		3586P	10 cm x 10 cm	5 dressings
3589P		10 cm x 15 cm	5 dressings	
3M Tegagen™ Alginate Dressing	90110	5 cm x 5 cm	10 dressings	
	90112	10 cm x 10 cm	10 dressings	
	90114	10 cm x 20 cm	5 dressings	
	90120	30,4 cm x 2 cm	5 dressings	
3M Tegaderm™ Transparent Dressing	1622W	4,4 cm x 4,4 cm	100 dressings	
	1624W	6 cm x 7 cm	100 dressings	
	1626W	10 cm x 12 cm	50 dressings	
	1627	10 cm x 25 cm	20 dressings	
	1628	15 cm x 20 cm	10 dressings	
	1629	20 cm x 30 cm	10 dressings	

Note about disposal

These tapes and dressings can be disposed of in the same way as household waste.



The skin is a complex and vital organ, which must be kept healthy; not only free from disease or injury, but also in a state of optimal well-being. 3M Skin Health is the specific area of 3M Health Care focused on both chronic and acute care. 3M Skin Health is committed to maintaining healthy skin, preventing the risk of skin disease and treating skin disease.



3M Health Care Limited

3M House, Morley Street
Loughborough
Leicestershire LE11 1EP
Tel: (01509) 613467
Fax: (01509) 613194
www.3M.com/healthcare/UK

3M Ireland

Adelphi Centre
Dun Laoghaire
Co Dublin, Ireland
Tel: (01) 2803 555
Fax: (01) 2803 509

3M, Durapore, Transpore, Micropore, Blenderm, Microfoam, Medipore, Tegaderm, Tegagen are trademarks of the 3M Company.

EF 519 0903 09 Graphix/VA/0903



The complete range ...

3M Health Care

innovation in action

3M™ Surgical Tapes & Dressings

Cover, care and protect


... provide the right surgical tape and adhesive dressing for every situation.

The skin is a delicate and sensitive organ. Patients with skin problems are especially at risk. So even small things – such as take right choice of surgical tapes and dressings – quickly take on greater importance.

3M is the world's largest manufacturer of surgical tapes and transparent dressings. 3M products are the result of the most up-to-date international research – they protect and care for the skin.

- Porous and breathable backing materials ensure the skin can breathe normally.
- All 3M surgical tapes and dressings are latex free.
- All 3M surgical tapes and dressings are pressure sensitive.
- Skin-friendly hypoallergenic polyacrylate adhesive reduces skin irritation.
- Selective adhesive application enables optimal moisture exchange while providing excellent adhesion.

Overview of Average Adhesion of 3M Surgical Tapes & Dressings

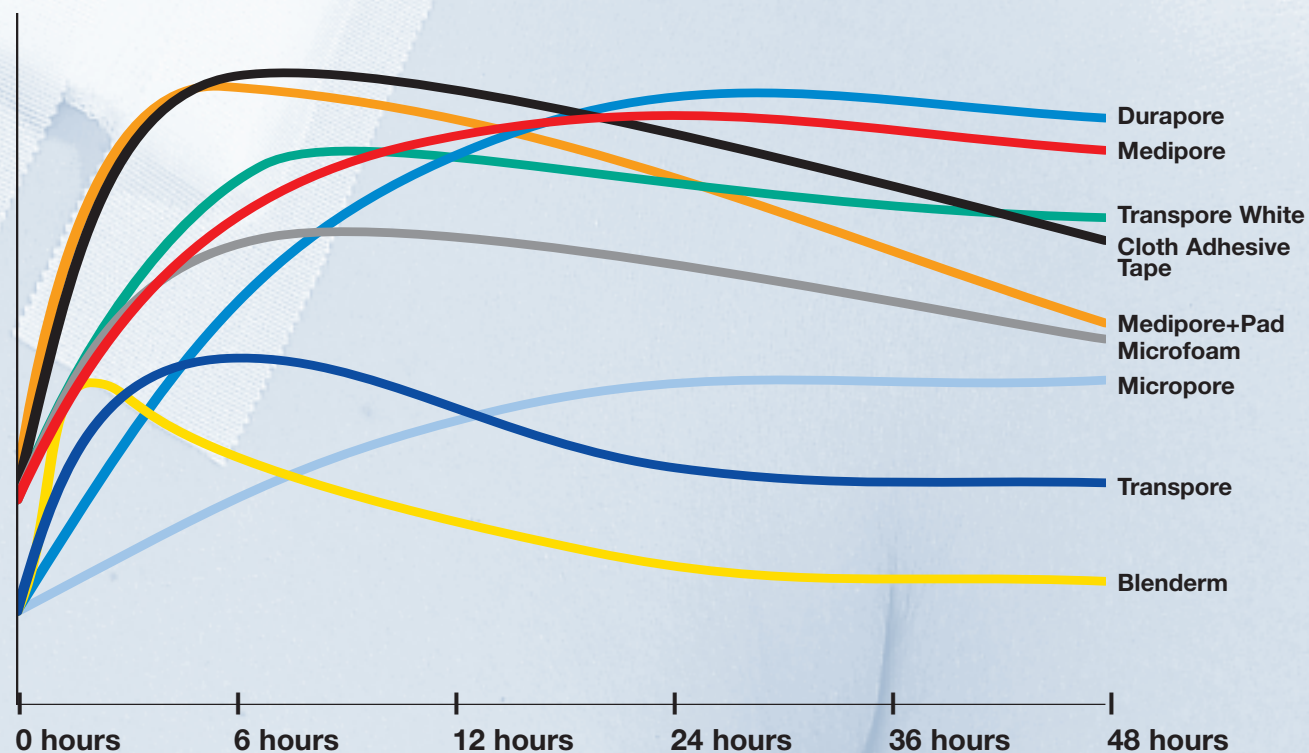


Successful Tape Application



Apply the tape evenly by positioning and pressing down from the centre towards the edges.



Optimise adhesion by firmly rubbing or pressing the tape backing. This will ensure the tape stays in place.

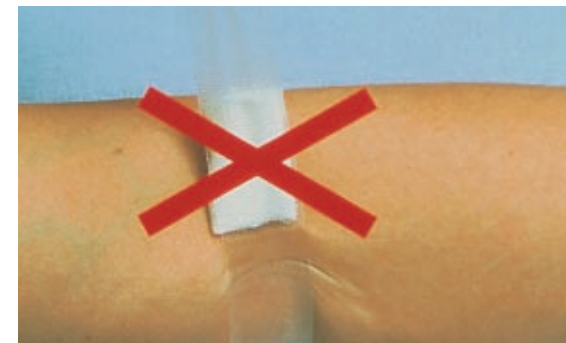

Never apply the tape under tension. This may induce skin trauma in the form of blistering or skin tearing.

Successful Removal



When removing, loosen ends and remove in the direction of hair growth.




Keep the tape close to the skin surface as you slowly pull it back and support the skin throughout removal.

Good Taping: Choosing the Right Tape or Dressing

Applications

3M products recommended

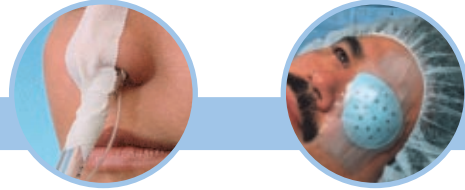

General secure when stretch is not needed

- Securing tubing (i.e. endotracheal and nasogastric tubes, chest tubes, I.V.)
3M Cloth Adhesive Tape, 3M Durapore™ Tape, 3M Transpore™ Tape, 3M Transpore™ White Tape
- Securing catheters e.g. urinary, newly inserted central I.V. catheters
3M Transpore Tape, 3M Medipore™ + Pad, 3M Tegaderm™ + Pad
- Securing other devices e.g. colostomy bags
3M Durapore Tape, 3M Micropore™ Tape, 3M Medipore™ H Tape, 3M Medipore™ Tape
- Adhering and covering simple dressings
3M Cloth Adhesive Tape, 3M Micropore Tape, 3M Transpore White Tape, 3M Transpore Tape
- Anchoring bulky dressings e.g. eye dressings
3M Transpore Tape, 3M Durapore Tape



Immobilization

- Taping bandages (strapping)
3M Cloth Adhesive Tape
- Immobilizing fingers and toes
3M Cloth Adhesive Tape, 3M Durapore Tape, 3M Transpore White
- Splints fixation
3M Cloth Adhesive Tape, 3M Durapore Tape




High Risk or Delicate Skin

- Taping applications on fragile skin including facial dressings, repeated taping applications, eye covers/dressings
3M Micropore, 3M Transpore White
- When repeated taping is needed where preserving "at-risk" skin is critical i.e. chronic wound dressings
3M Medipore H Tape, 3M Medipore Tape with liner

- **The more critical for tubing/dressing, the higher the adhesion required**
- **Strength of adhesion and backing may be more important than gentleness**
- **Gentle tapes with higher adhesion may be used without increasing the risk of skin trauma if properly applied and removed**
- **Some tapes are less gentle (e.g. Durapore tape and Cloth Adhesive tape) but are used for high initial adhesion, high long-term adhesion, high strength of backing.**

Applications

3M products recommended

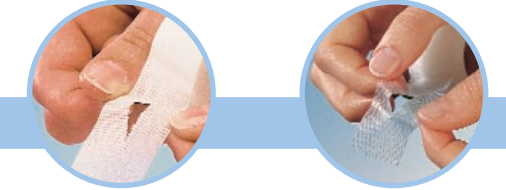

When tape width must be customized

3M Transpore White Tape, 3M Transpore Tape



When you need a waterproof tape or dressing

- Protecting dressings from external moisture, fluids or contaminants
3M Blenderm™ Tape
- Treating dermatological disorders
- Adjunct to the treatment of dermatological conditions
- Occlusive test patches
- Protecting newly inserted central I.V. catheter sites
3M Tegaderm + Pad



When you need a tape that stretches

(e.g. anticipate edema, distention, hematoma formation or movement)

- Adhering and covering post-surgical dressings: abdominal, cardio-thoracic and orthopaedic, OB/GYN
- For compression over cardiovascular/angio sites, breast surgery and other surgical sites
3M Medipore H Tape, 3M Microfoam™ Tape



When you need an absorbent dressing

- Primary or secondary dressing over light to moderately draining wounds including: surgical wounds, traumatic wounds such as cuts, burns, abrasions
3M Medipore + Pad, 3M Tegaderm + Pad
- Secondary dressing over hydrogels, alginates and gauze
3M Tegaderm + Pad
- Primary dressing over high exuding wounds including: total hip implant surgery wound (based on clinical evidence*), thoracic wounds
- Reduction of skin blistering in cases of total hip implant surgery
- Prolonged wear time (up to 12 days*)
*study performed on 299 patients in Norway

Overview

3M Surgical Tapes and Dressings



	Cloth Adhesive Tape	Durapore Tape	Transpore White Tape	Transpore Tape	Micropore Tape	Blenderm Tape
	Hypoallergenic, strong adhesion	Silk-like cloth tape, universal	Gentle tape with strong adhesion, easy-to-tear	Transparent plastic tape, easy-to-tear	Skin-friendly	Waterproof
Backing material	Cotton cloth	Viscose Silk taffeta	Nonwoven	Ethylene vinyl acetate	Viscose, nonwoven	Polyethylene film
Adhesive	Acrylate Adhesive Zinkoxyde	Polyacrylate	Polyacrylate	Polyacrylate	Polyacrylate	Polyacrylate
Hypoallergenic	X	X	X	X	X	X
Porous		X			X	
Breathable	X	X	X	X	X	
Strong backing	X	X				
Easy to use with gloves	X					
Easy to tear off	X		X	X	X	
Bi-directional tear			X	X		
Minimal adhesive residue upon removal from skin	X	X	X	X	X	X
Initial adhesion	High on damp and dry skin		Holds well on damp and dry skin	Excellent		
Can be written on	X					
Transparent				X		X
Waterproof						X
Conformable			X			X
Sticks well to itself			X			
Soft					X	
Occlusive						X
Perforated						
Liner						
Sterile						
Absorbent pad						
"Frame-delivery" system						
Low adherent						
Bacterial and viral barrier						



	Microfoam Tape	Medipore H Tape Perforated	Medipore Tape on liner	Medipore + Pad Dressing	Tegaderm + Pad Dressing	Tegaderm Dressing	Tegagen Alginate Dressing
	Elastic foam tape, stretchy	Conformable tape with "easy-tear" perforations	Conformable tape with liner	Quickly and easily applied	Shower- proof dressing	Transparent wound dressing	Highly absorbent dressing
Backing material	Polyvinyl chloride foam, nonwoven	Polyester nonwoven	Polyester nonwoven	Polyester nonwoven	Polyurethane viscose	Polyurethane film	high molecular weight hydrophilic polymer which is extracted from sea weed.
Adhesive	Polyacrylate	Polyacrylate	Polyacrylate	Polyacrylate	Polyacrylate	Polyacrylate	-
Hypoallergenic	X	X	X	X	X	X	
Porous		X	X	X			X
Breathable	X	X	X	X	X	X	X
Strong backing							
Easy to use with gloves							
Easy to tear off		X					
Bi-directional tear							
Minimal adhesive residue upon removal from skin							
Initial adhesion							
Can be written on							
Transparent					X	X	
Waterproof	X				X	X	
Conformable	X	X	X	X	X	X	X
Sticks well to itself							
Soft		X	X	X			
Occlusive							
Perforated		X					
Liner			X				
Sterile				X	X	X	X
Absorbent pad				X	X		
"Frame-delivery" system					X	X	
Low adherent				X	X		
Bacterial and viral barrier					X	X	